


# APSW 4-PLY, WIRE-REINFORCED, MANDREL-WRAPPED, SILICONE SUCTION HOSE

PLATINUM-CURED; SMOOTH OR CONVOLUTED O.D.

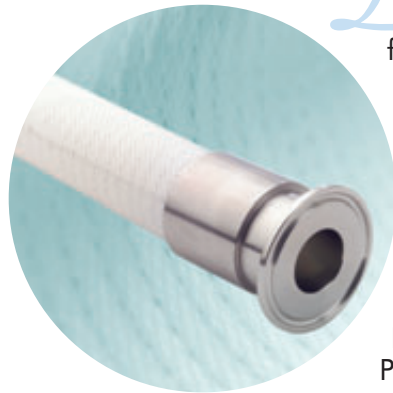

Low volatile grade, platinum-cured silicone hose, produced in our certified **Class 1000 clean room** for critical pharmaceutical, biomedical, cosmetic, and food applications, comprises the core of APSW. The core is then slipped onto a mandrel and wrapped with polyester mesh fabric, 316 stainless steel reinforcing wire, and additional silicone. It is cured into a homogeneous hose to enhance its pressure and vacuum capabilities. APSW has undergone extensive physical, chemical, and biological testing and meets USP **Class VI**, **FDA CFR 177.2600**, ISO 10993, European Pharmacopoeia 3.1.9, and 3-A standards.

APSW is not intended for implantation.



## KEY FEATURES



Core


- Available with either a smooth or convoluted O.D.
- Low volatile grade silicone suitable for pharmaceutical, biomedical, cosmetic, and food applications such as ultrapure water transfer, food processing, washdown stations, skid systems, and other suction and discharge uses
- Core made of **NSF 51** listed silicone (National Sanitation Foundation)
- Hardness value of 70 Shore A
- Rated for full vacuum
- Crush and kink resistant
- Durable and heavy-duty yet translucent for visual flow contact
- Two styles available: Smooth O.D. or, for an improved range of flexibility, Convoluted
- General temperature range: -100°F (-73.3°C) to 400°F (204.4°C)
- Sterilizable by autoclave, CIP, SIP, and gamma radiation processes
- Resists temperature extremes, chemical attack, ozone, radiation, moisture, and environmental exposure
- Imparts no taste or odor to critical streams
- Stocked in 12 ft. straight lengths
- Manufactured under strict Good Manufacturing Procedures in a tightly-controlled environment
- Wrapped in clear polyethylene at manufacture to ensure cleanliness
- Documented lot traceable with identification on bags and boxes
- Documented quality control
- Complete validation package available upon request
- Custom lengths, sizes, special cleaning, and/or packaging available
- Available with Nomex® braid construction – call for details

SILICONE

[www.advantapure.com](http://www.advantapure.com)

# APSW 4-PLY, WIRE-REINFORCED, MANDREL-WRAPPED, SILICONE SUCTION HOSE


PLATINUM-CURED; SMOOTH OR CONVOLUTED O.D.

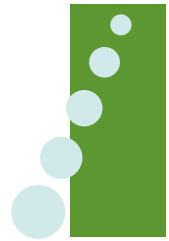

## SPECIFICATIONS


Product Number	I.D.		O.D.		Vacuum Rating (Hg.)	Working Pressure	Burst Pressure*		Min. Bend Radius	Weight per	
	(in.)	(mm)	(in.)	(mm)	at 130°F (54.4°C)	at 70°F (21.1°C) (PSI) (Bar)	at 70°F (21.1°C) (PSI) (Bar)	(in.)	Foot / Meter	(lb)	(kg)

### APSW-P SMOOTH O.D.

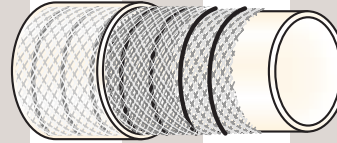

APSW-P-0500	.500	12.70	.910	23.11	29.9	150	10.3	600	41.4	2.0	.281	.418
APSW-P-0750	.750	19.05	1.16	29.46	29.9	150	10.3	600	41.4	2.5	.359	.534
APSW-P-1000	1.000	25.40	1.41	35.81	29.9	150	10.3	600	41.4	3.5	.438	.652
APSW-P-1500	1.500	38.10	1.91	48.51	29.9	150	10.3	600	41.4	4.0	.685	1.019
APSW-P-2000	2.000	50.80	2.41	61.21	29.9	125	8.6	500	34.5	6.0	.875	1.302
APSW-P-2500	2.500	63.50	2.91	73.91	29.9	125	8.6	500	34.5	11.0	1.050	1.562
APSW-P-3000	3.000	76.20	3.41	86.61	CF	100	6.9	400	27.6	13.0	1.230	1.830
APSW-P-4000	4.000	101.60	4.41	112.01	CF	100	6.9	400	27.6	CF	CF	CF

### APSW-PC CONVOLUTED O.D.

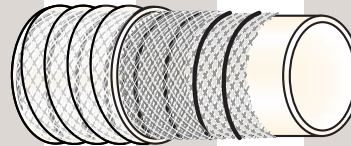

APSW-PC-0500	.500	12.70	.910	23.11	29.9	150	10.3	600	41.4	2.0	.208	.309
APSW-PC-0750	.750	19.05	1.16	29.46	29.9	150	10.3	600	41.4	2.5	.307	.457
APSW-PC-1000	1.000	25.40	1.41	35.81	29.9	150	10.3	600	41.4	3.5	.390	.580
APSW-PC-1500	1.500	38.10	1.91	48.51	29.9	150	10.3	600	41.4	4.0	.625	.930
APSW-PC-2000	2.000	50.80	2.41	61.21	29.9	125	8.6	500	34.5	6.0	.775	1.153
APSW-PC-2500	2.500	63.50	2.91	73.91	29.9	125	8.6	500	34.5	CF	CF	CF
APSW-PC-3000	3.000	76.20	3.41	86.61	CF	100	6.9	400	27.6	CF	CF	CF
APSW-PC-4000	4.000	101.60	4.41	112.01	CF	100	6.9	400	27.6	CF	CF	CF

Wall dimensions for all sizes and styles: .205"; 5.21mm

CF = Consult Factory

\*For every 100°F of temperature over 70°F, reduce the burst pressure by 5%.

**NOTE:** When products are used as part of an assembly, the pressure ratings of fittings may be less than hose pressure ratings above. Please consult your sales representative.

SILICONE

Purity in Fluid Flow Systems<sup>SM</sup>



Address: NewAge® Industries AdvantaPure®  
145 James Way  
Southampton, PA 18966 U.S.A.

**AdvantaPure**  
NEWAGE INDUSTRIES



DISTRIBUTED BY:

Toll Free Phone: **888-755-4370**  
Phone: **215-526-2151**  
Toll Free Fax: **888-258-4293**  
Fax: **215-526-2167**  
Web Site: **www.advantapure.com**  
E-mail: **sales@advantapure.com**

©2006 NewAge® Industries, Inc.  
AdvantaPure products are produced by NewAge® Industries, Inc.  
NewAge® Industries AdvantaPure®, the AdvantaPure logo, NewAge® and the "N" logo are registered trademarks of NewAge® Industries, Inc.  
Purity in Fluid Flow Systems<sup>SM</sup> is a service mark of NewAge® Industries, Inc.  
Nomex® reg TM E. I. DuPont de Nemours and Company

WARRANTY: NewAge Industries AdvantaPure products are warranted to be free from defects in material and workmanship. Any product found to carry such defects will be replaced without charge. This warranty is limited only to the replacement of NewAge Industries AdvantaPure products, item for item, which after examination by NewAge Industries AdvantaPure are deemed defective. NewAge Industries AdvantaPure specifically disclaims any other liability. The data supplied is provided as a helpful guideline and is believed to be reliable; however, nothing stated shall constitute a guarantee, recommendation, or warranty for any application.

NOTICE: NewAge Industries, Inc. has been a privately owned corporation since 1954. Neither the corporation nor its principals are affiliated, in any way, with any other group or entity either representing themselves as or being referred to by the term New-Age.



**ESOP** Employee Owned  
for Your Benefit

Celebrating over  
**50**  
Years  
of Service  
1954-2006